

## Neck for Tracheal Stenosis 3.5

**Indication:** Tracheal Stenosis

**Patient Preparation:** Remove metallic dentures and jewelry

**IV Contrast:** None


Start at Sub-Mandibular through the carina

Scan Mode	Thickness Speed Pitch	Table Interval	SFOV	kVp	Auto mA/ Noise Index	Prep Time (sec)	Recon Type
VCT 0.8 sec Helical Full	1.25 9.37 0.984:1	1.25	Large 23	120	440 6	N	Bone Standard
LS 16 0.8 sec Helical Full	1.25 9.37 0.984:1	1.25	Large 23	120	440 6	N	Bone Standard
AS 64 1	1	1	23	120	300	N	Bone Standard

**Algorithm:** Standard and bone. Recon (1) 1.25 Standard (2) 1.25 Bone (3) 0.625 Standard

**Network:** Recon 1 and 2 to PACS. Recon 1.2 to AWMR

**Comments:** Coronal and Sagittal reformats; 3D Volume Rendering. (Need to order Image Processing with Workstation).

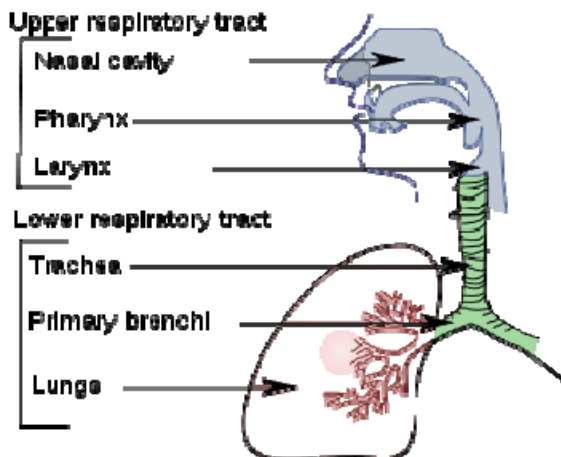

Written: 07/2010

LSUHSC Shreveport  
Radiology Department  
Proc 13.3.23